

Efficient Ion Fragmentation in Structures for Lossless Ion Manipulations



Pacific Northwest
NATIONAL LABORATORY

Proudly Operated by Battelle Since 1965

Ian K. Webb, Sandilya Garimella, Isaac K. Attah, Christopher D. Chouinard, Gabe Nagy, Aneesh Prabhakaran, Erin S. Baker, Yehia M. Ibrahim and Richard D. Smith
Pacific Northwest National Laboratory, Richland, WA

Overview

- Fragmentation of peptide cations has been achieved in structures for lossless ion manipulations (SLIM).

- A DC offset was applied to RF electrodes on one of two parallel SLIM surfaces.

- b and y-type sequence fragment ions were produced.

Introduction

- SLIM have been used for ultra high resolution ion mobility separations.¹

- Flexibility of design and continuous RF confinement allow complex 3D ion manipulations over extended times.²

- SLIM are being applied to separations and quantitation of complex omics mixtures.

- Identification of unknown compounds without internal standards is difficult.

- Mass spectrometric fragmentation (e.g. collision induced dissociation, CID) allows for identification of unknown compounds.

- Therefore, ion fragmentation in SLIM has been developed.

- Ion trajectories were simulated and electric fields were modeled to assist in the inference of the activation mechanism.

Methods

- Micromolar concentration peptides and proteins were infused by nano-electrospray ionization at 300 nl/min at +3000 V.

- DC bias was applied to the RF of the top SLIM surface via a 100 kΩ coupling resistor to fragment ions.

- 30 V, 200 m/s traveling wave with 1 MHz, 250 V_{pp} RF

- SLIM was held at 3 Torr N₂.

- Simulations were performed in Simion 8.1 with the SDS collision model.


Results

Simulated Ion Trajectories and Kinetic Energies for Triply Protonated Neurotensin




0 V Bias


25 V Bias

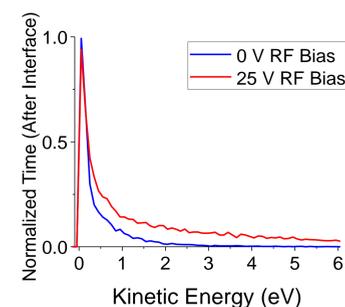


CID of Neurotensin Peptide Cations

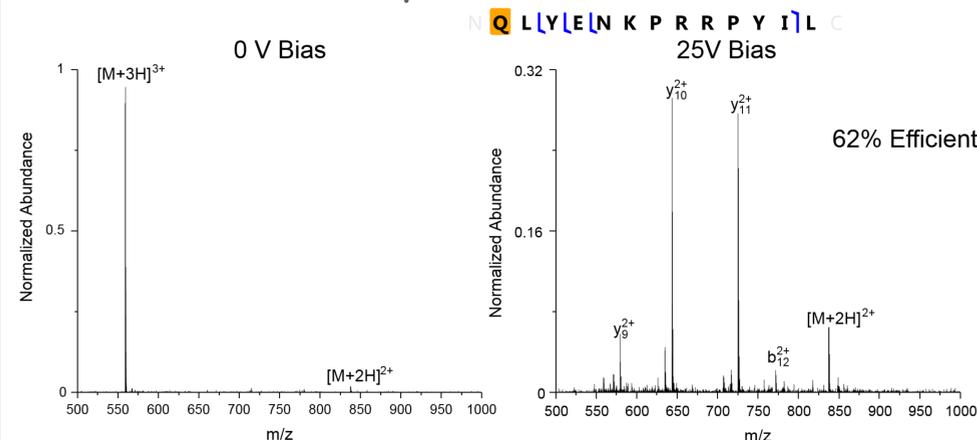

CID of Melittin Peptide Cations

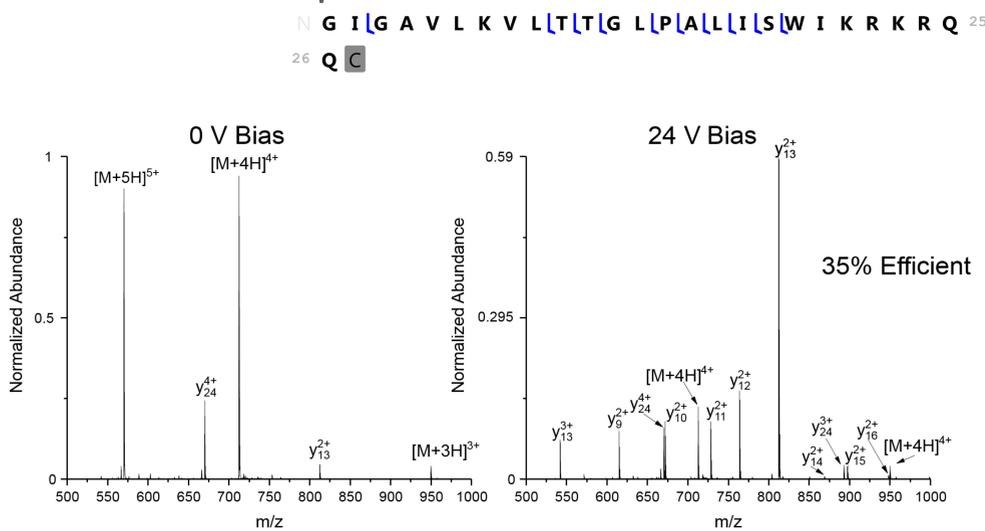

CID of Ubiquitin Cations

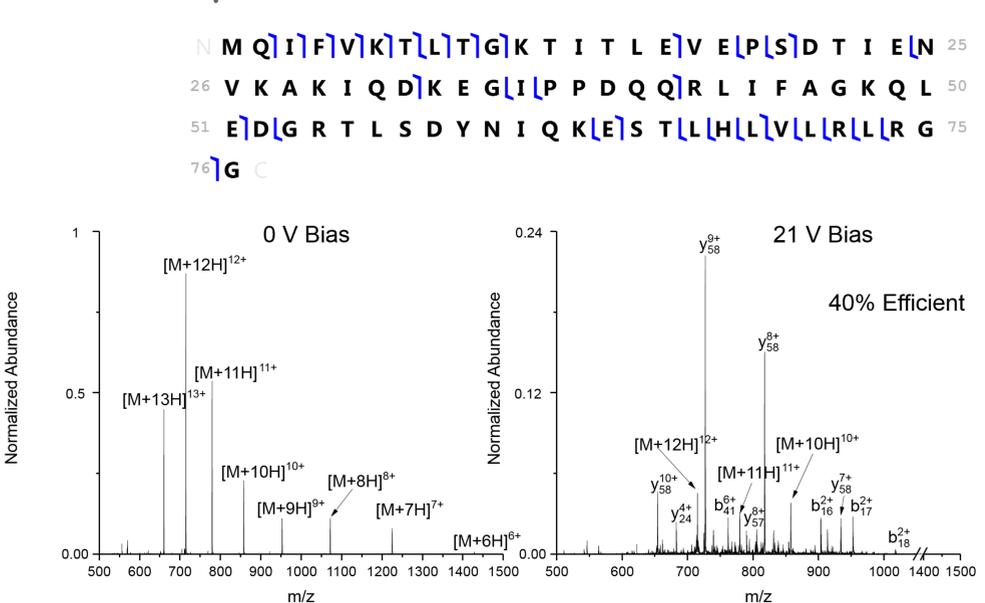

Conclusions

- A simple method for implementing CID with SLIM has been demonstrated.

- Ion activation was induced by applying a DC bias to the RF on a single SLIM surface..

- Simulations of ion trajectories show ions traveling into regions of higher RF field intensities when bias is applied.

- Resulting ion kinetic energies are higher when bias is applied.

- We plan to implement CID with high resolution ion mobility separations, accumulation of large ion populations, and other manipulations.

Acknowledgements

This work was supported by the NIH National Institute of General Medical Sciences (8 P41 GM103493-10), the U.S. Department of Energy (DOE) Office of Biological and Environmental Research (OBER) Pan-omics program at Pacific Northwest National Laboratory (PNNL), and MOBILion Systems Inc. Work was performed in the Environmental Molecular Sciences Laboratory, a DOE OBER national scientific user facility on the PNNL campus. PNNL is a multiprogram national laboratory operated by Battelle for the DOE under contract DE-AC05-76RL01830.

References

- Deng, L., Webb, I.K., et al.: Serpentine ultralong path with extended routing (SUPER) high resolution traveling wave ion mobility-MS using structures for lossless ion manipulations. *Anal. Chem.* **89**(8), 4628-4634, (2017)
- Ibrahim, Y.M., et al.: New frontiers for mass spectrometry based upon structures for lossless ion manipulations. *Analyst*, **142**(7), 1010-1021, (2017)

CONTACT: Ian K. Webb, Ph.D.
Biological Sciences Division
Pacific Northwest National Laboratory
E-mail: ian.webb@pnnl.gov

Career Opportunities: For potential openings in Integrative Omics Group at PNNL please visit <http://omics.pnl.gov/careers>

www.omics.pnl.gov

